


# CUTIMED® DRESSING SELECTION GUIDE

THIS CUTIMED® DRESSING SELECTION GUIDE HELPS YOU TO CHOOSE THE RIGHT DRESSING FOR INFECTED, EXUDING WOUNDS.

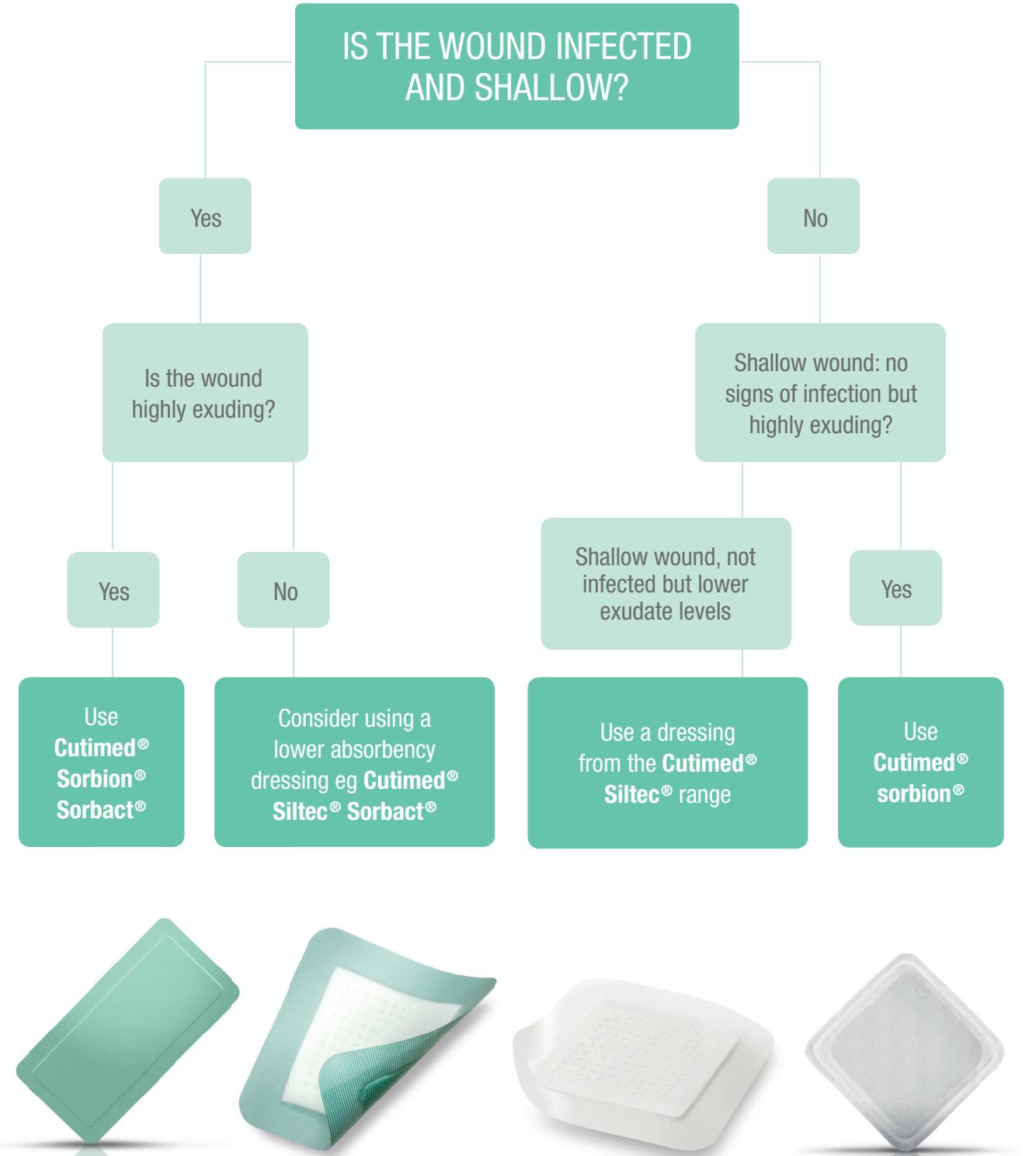

## EFFICACY BASED ON EVIDENCE AND EXPERIENCE

Cutimed® wound care products have been proven safe and effective in the management of difficult-to-heal wounds. Please visit: [www.cutimed.com/awc/evidence/](http://www.cutimed.com/awc/evidence/) to see the complete list of downloadable publications and hands-on case reports.

In addition, Cutimed® wound care products are also effective in therapy approaches specifically developed for venous leg ulcers, diabetic ulcers and pressure ulcers. Each of these therapy options can be downloaded at: <http://cutimed.com/therapy-approaches/> or by contacting BSN medical.



BSN medical GmbH

Quickbornstraße 24  
20253 Hamburg

[www.cutimed.com](http://www.cutimed.com)  
[www.bsnmedical.com](http://www.bsnmedical.com)  
Tel +49 40 / 49 09-909  
Fax +49 40 / 49 09-6666

70551-000008-00 07/2015

**BSN** medical

THERAPIES. HAND IN HAND.

**Cutimed® Sorbion® Sorbact®**

Closing wounds. Together.

DRAFT



**BSN** medical

# CUTIMED® SORBION® SORBACT® A UNIQUE COMBINATION

All wounds can present challenges to healthcare professionals, but infection and very high levels of exudate can further complicate their management. A dressing is required that is both antimicrobial and absorbent and which can safely lock away wound fluid. Treatment options for infected wounds often include silver based dressings, but these can cause side effects or may not be suitable for all patient groups.

To help care providers with these challenging wounds, BSN medical has developed the first antimicrobial, super-absorbent wound dressing: **Cutimed® Sorbion® Sorbact®**.

Combining the proven antimicrobial properties of **Cutimed® Sorbact®** with **super-absorbent technology**, this unique dressing provides reliable exudate management plus bioburden reduction without any contraindications.



## Indications

**Cutimed® Sorbion® Sorbact®** is recommended for all shallow, clean, colonised, contaminated or infected and highly exuding wounds including:

- Venous leg ulcers,
- Diabetic foot ulcers,
- Dehisced surgical wounds.


It can be also used as a secondary absorbent dressing on deeper wounds.

# ONE DRESSING FOR INFECTED AND HIGHLY EXUDING WOUNDS

Each layer of **Cutimed® Sorbion® Sorbact®** is designed to provide an optimal treatment outcome:

## HYDROPHOBIC NON-WOVEN BACKING


- Low risk of strike-through
  - ▶ Protects patients clothes and bedding
  - ▶ Improving patient comfort
- Protects from external contamination
  - ▶ Wound healing is undisturbed
- Printed surface
  - ▶ Enables easy dressing application

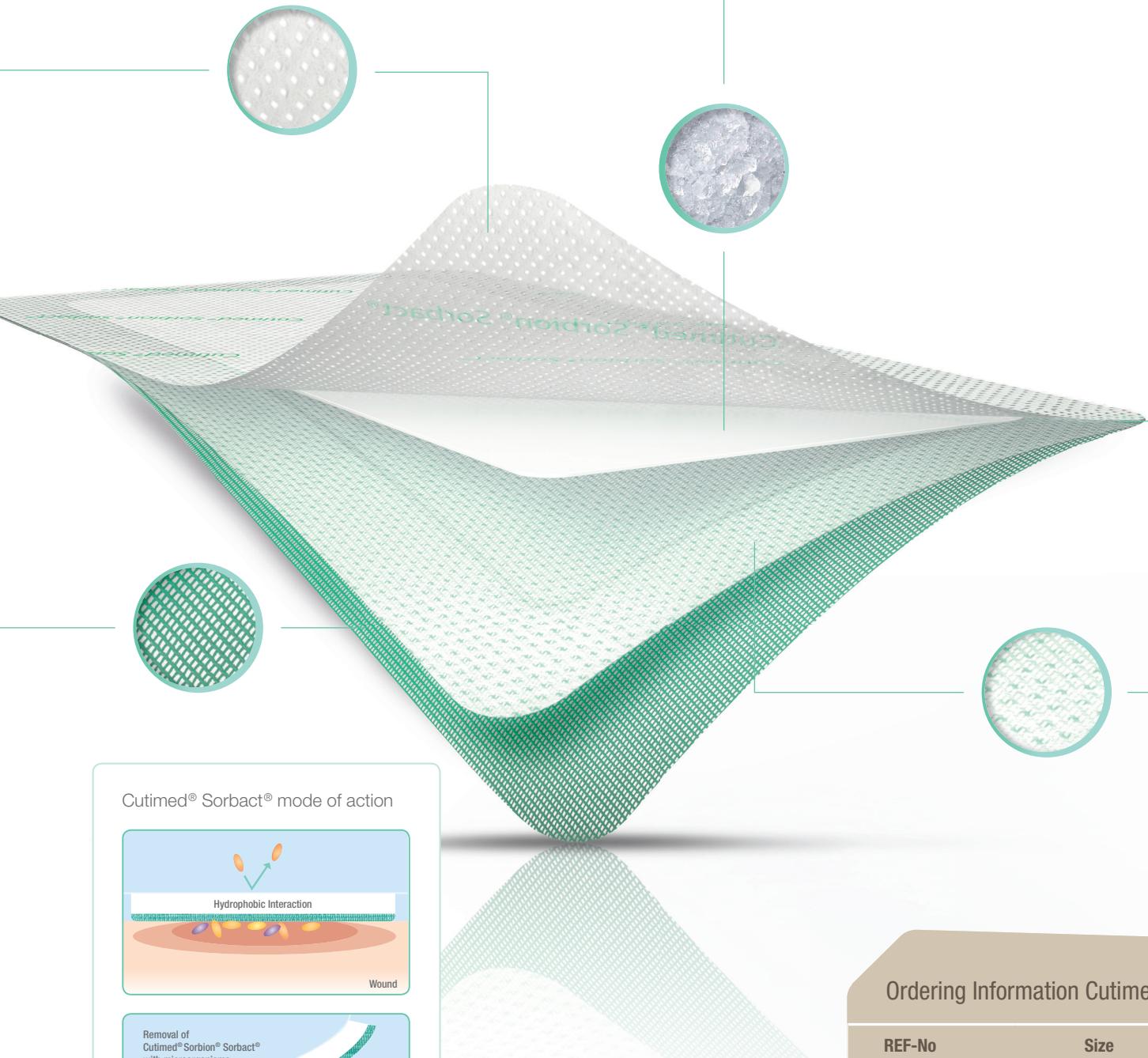

## CUTIMED® SORBACT® WOUND CONTACT LAYER

Due to its coating of DACC, **Cutimed® Sorbact®** enables safe<sup>1</sup>, irreversible, physical binding of bacteria and fungi to the dressing and rendering them inert<sup>2</sup>.

- Low risk of allergies or contraindications
  - ▶ Can be used safely on all patients including babies, pregnant women, nursing mothers and those previously sensitised to other antimicrobial dressings
- No use of chemically active agents
  - ▶ No bacterial or fungal resistance or contraindications
- Leaves no cell debris behind
  - ▶ Supporting natural wound healing

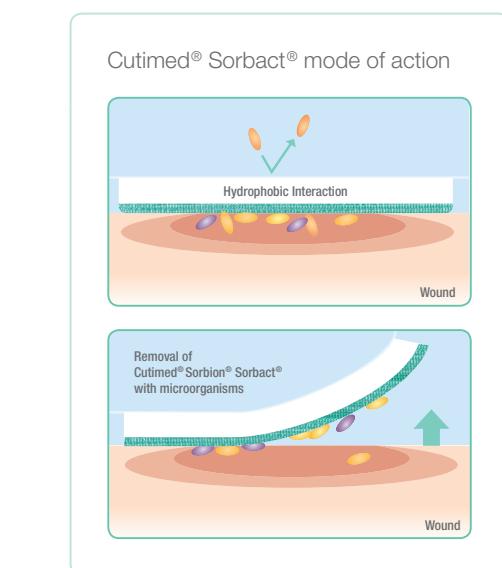


<sup>1</sup> Haycocks S, Chadwick P (2011). Use of a DACC coated antimicrobial dressing in people with diabetes and a history of foot ulceration. Wounds UK Vol 6 No 4

<sup>2</sup> Lungh et al (2006) Using the principle of hydrophobic interaction to bind and remove wound bacteria. Journal of Wound Care, 15 (4): 175-80

## SUPER-ABSORBENT CORE

- Absorbs and retains large volumes of exudate into the dressing even under pressure
- ▶ Prevents skin maceration and assists with the management of different wounds e.g. venous leg ulcers and pressure ulcers

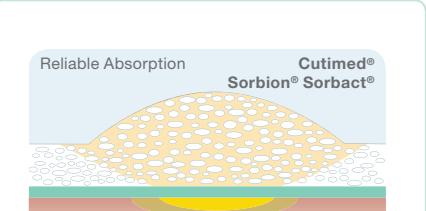

The **Cutimed® Sorbion® Sorbact®** internal core expands and gels to effectively retain high to very high levels of exudate



Cushioning gelling core protects the wound and offers increased patient comfort

## ROUNDED EDGES

- Remain flat
- ▶ Provides additional patient comfort

## NON-WOVEN LAYER

- Allows for optimal distribution of fluid throughout the dressing and prevents exudate returning to the wound bed
- ▶ Reduces risk of skin maceration

## Ordering Information Cutimed® Sorbion® Sorbact®

REF-No	Size	Wound pad size	Box of	Shipping Unit
72698-00	10 cm x 10 cm	8 x 9 cm	10	12 boxes
72698-01	10 cm x 20 cm	7.7 x 17.8 cm	10	6 boxes
72698-02	20 cm x 20 cm	17.8 x 17.8 cm	10	4 boxes
72698-03	20 cm x 30 cm	17.5 x 27.5 cm	10	4 boxes